

## DEVELOPMENT OF PCR TEST FOR DETECTION OF THE SARS-COV-2 GENETIC VARIANTS ALPHA, BETA, GAMMA, DELTA

Shipulin GA, Savochkina YuA, Shuryaeva AK , Shivyagina EE, Nosova AO, Davydova EE, Luparev AR, Malova TV, Yudin SM


Centre for Strategic Planning and Management of Biomedical Health Risks of Federal Medical Biological Agency, Moscow, Russia

The emergence of novel SARS-CoV-2 genetic variants with increased transmissivity and reduced antibody neutralization efficiency is a threat to global public health. Reverse transcription polymerase chain reaction (RT-PCR) with the use of fluorescent probes, which make it possible to detect the single nucleotide substitutions, is a technique suitable for screening the SARS-CoV-2 RNA-containing samples for the already known functionally significant mutations in the S-gene, identification of which allows to define and differentiate the most epidemiologically significant genetic variants. The study was aimed to develop an assay for the large-scale monitoring of the spread of the SARS-CoV-2 top-priority variants. Based on the whole-genome alignment of the SARS-CoV-2 sequences, deposited in the GISAID database, primers and LNA-modified probes were selected to detect mutations in the S gene, typical for the Alpha, Beta/Gamma and Delta variants of concern (VOC). The developed reagent kit for detection of the key mutations in the SARS-CoV-2 S gene by the real time RT-PCR has good analytical and diagnostic characteristics and was authorized as a medical device (reagent) for *in vitro* use. The results of detecting the VOC and the key mutations with the use of the developed reagent kit were consistent with the data of the whole genome sequencing of 1,500 SARS-CoV-2 RNA samples. The developed reagent kit and the subsequent SARS-CoV-2 RNA sequencing assay used to perform the epidemiological monitoring of SARS-CoV-2 variants made it possible to promptly report the emergence of the Delta genetic variant in Russia, and to trace the dynamic changes in the prevalence of Delta in Moscow Region in April–September 2021.

**Keywords:** coronavirus, COVID-19, SARS-CoV-2, N501Y, P681H, 69-70del, E484K, B.1.1.7, B.1.351, P.1, VOC

**Author contribution:** Shipulin GA, Savochkina YuA, Davydova EE, Yudin SM — planning the experiment; Shipulin GA, Savochkina YuA, Yudin SM — literature analysis; Shuryaeva AK, Shivyagina EE, Nosova AO, Luparev AR, Malova TV — experimental procedure, data interpretation; Savochkina YuA, Shivyagina EE, Nosova AO — reagent kit development; Luparev AR — statistical analysis; Shipulin GA, Shuryaeva AK — manuscript writing and editing; Davydova EE, Yudin SM — manuscript editing.

**Compliance with ethical standards:** the study was performed in accordance with the requirements of the Declaration of Helsinki and GOST R ISO 14155-2014.

 **Correspondence should be addressed:** Anna K. Shuryaeva  
Pogodinskaya, 10, str. 1, 119121, Moscow; ashuryaeva@cspmz.ru

**Received:** 15.12.2021 **Accepted:** 13.01.2022 **Published online:** 09.02.2022

**DOI:** 10.47183/mes.2022.003

## РАЗРАБОТКА ПЦР-ТЕСТА ДЛЯ ВЫЯВЛЕНИЯ ГЕНЕТИЧЕСКИХ ВАРИАНТОВ АЛЬФА, БЕТА, ГАММА, ДЕЛЬТА ВИРУСА SARS-COV-2

Г. А. Шипулин, Ю. А. Савочкина, А. К. Шуряева , Е. Е. Шивлягина, А. О. Носова, Е. Е. Давыдова, А. Р. Лупарев, Т. В. Малова, С. М. Юдин


Центр стратегического планирования и управления медико-биологическими рисками здоровья Федерального медико-биологического агентства, г. Москва, Россия

Возникновение новых вариантов коронавируса SARS-CoV-2, обладающих повышенной трансмиссивностью и снижающих эффективность его нейтрализации выработанными ранее антителами, представляет угрозу для здоровья населения во всем мире. Метод полимеразной цепной реакции с обратной транскрипцией (ОТ-ПЦР) с флуоресцентными зондами, позволяющими детектировать единичные нуклеотидные замены, подходит для скрининга содержащих РНК SARS-CoV-2 образцов на наличие уже известных мутаций в S-гене, которые имеют функциональное значение, и выявление которых позволяет определять и дифференцировать геноварианты, имеющие наибольшее эпидемиологическое значение. Целью настоящей работы было разработать набор реагентов и методику для оперативного мониторинга распространения приоритетно значимых вариантов SARS-CoV-2. На основе выравнивания полногеномных последовательностей SARS-CoV-2, опубликованных в базе данных GISAID, были выбраны праймеры и LNA-модифицированные зонды для выявления мутаций в S-гене, характерных для генетических линий альфа, бета/гамма и дельта особой эпидемиологической значимости (VOC). Разработан и зарегистрирован набор реагентов в формате ОТ-ПЦР в реальном времени для выявления наиболее значимых мутаций в S-гене SARS-CoV-2, продемонстрированы его высокие аналитические и диагностические характеристики. Показано соответствие результатов выявления линий VOC и их основных мутаций разработанным набором с данными полногеномного секвенирования для 1500 образцов РНК SARS-CoV-2. Применение набора реагентов в сочетании с последующим секвенированием РНК SARS-CoV-2 в рамках проведения эпидемиологического мониторинга позволило оперативно установить факт появления на территории России геноварианта дельта и в дальнейшем отследить динамику изменения его распространенности в Московском регионе в период с апреля по сентябрь 2021 г.

**Ключевые слова:** коронавирус, COVID-19, SARS-CoV-2, N501Y, P681H, 69-70del, E484K, B.1.1.7, B.1.351, P.1, изоляты особой эпидемиологической значимости

**Вклад авторов:** Г. А. Шипулин, Ю. А. Савочкина, Е. Е. Давыдова, С. М. Юдин — планирование эксперимента; Г. А. Шипулин, Ю. А. Савочкина, С. М. Юдин — анализ литературы; А. К. Шуряева, Е. Е. Шивлягина, А. О. Носова, А. Р. Лупарев, Т. В. Малова — постановка эксперимента и интерпретация результатов; Ю. А. Савочкина, Е. Е. Шивлягина, А. О. Носова — разработка набора реагентов; А. Р. Лупарев — статистический анализ; Г. А. Шипулин, А. К. Шуряева — написание и редактирование рукописи; Е. Е. Давыдова, С. М. Юдин — редактирование рукописи.

**Соблюдение этических стандартов:** исследование проведено в соответствии с требованиями Хельсинкской декларации и ГОСТ Р ИСО 14155-2014.

 **Для корреспонденции:** Анна Константиновна Шуряева  
ул. Погодинская, д. 10, стр. 1, 119121, г. Москва; ashuryaeva@cspmz.ru

**Статья получена:** 15.12.2021 **Статья принята к печати:** 13.01.2022 **Опубликована онлайн:** 09.02.2022

**DOI:** 10.47183/mes.2022.003

A novel coronavirus, causing the dangerous respiratory disease in humans, COVID-19, was found at the end of the year 2019 in China. After the whole genome sequencing the virus was classified as a betacoronavirus and referred to as SARS-CoV-2 [1–3]. High mutagenic capability of coronaviruses is well known [4]. The majority of emerging mutations do not affect the properties of the virus, however, some mutations can result in functional alterations, including fast viral spread and/or more severe illness.

Since the start of SARS-CoV-2 pandemic, a large number of genome sequences were published in the GISAID database [5], which made it possible to trace both the viral spread and the emergence of mutations in the viral genome.

In December 2020, British authorities declared a rapid increase in the number of COVID-19 cases. The sharp spike in morbidity was caused by the novel SARS-CoV-2 coronavirus genetic variant, which differed from the reference genome. This variant carried mutations in the gene encoding the S protein [6, 7]. Based on the phylogenetic analysis, the new British lineage of SARS-CoV-2 was named B.1.1.7 and classified as a variant of concern (VOC). The B.1.1.7 variant has a number of typical mutations in the spike protein, such as N501Y, A570D, D614G, and P681H, as well as amino acid deletions 69-70 and 144Y [8, 9]. Simultaneously, one more epidemiologically significant lineage was detected, which caused the outbreak in various provinces of the Republic of South Africa. This lineage is known as B.1.351. The South African strain contains nine spike mutations in addition to D614G, including the mutation cluster (for example, 242–244del and R246I) in the spike protein N-terminal domain (NTD), three mutations (K417N, E484K and N501Y) in the receptor binding domain (RBD), and one mutation (A701V) located close to the furin cleavage site [10, 11]. In January 2021, experts of the National Institute of Infectious Diseases, Japan, discovered a novel SARS-CoV-2 variant, the P.1 lineage, the isolates of which were identified in tourists coming from Brazil. Some mutations present in the Brazilian strain were previously discovered in the British and South African variants. This isolate contains 12 mutations in the spike protein, including N501Y and E484K [12]. Later the discussed above coronavirus genetic lineages were named using the letters of the Greek alphabet in accordance with the nomenclature, introduced by the WHO: Alpha (B.1.1.7 and subvariants Q.1–Q.8), Beta (B.1.351), Gamma (P.1). All three SARS-CoV-2 genetic variants listed above are considered the variants of concern (VOC) by the WHO and ECDC, and circulate all over the world [6].

In the spring of 2021 in India the new SARS-CoV-2 lineage B.1.617.2 was discovered, which spread rapidly and gave rise to the new wave of the pandemic, with the sharp rise in the incidence rate and the number of deaths in India and then in other countries in different parts of the world. Along with some closely related genetic lineages identified most recently, AY.1 and AY.2-124, the B.1.617.2 lineage contains the following typical mutations in the S protein: T19R, E156G, del157/158, L452R, T478K, D614G, P681R, and D950N, including mutations in the RBD [13]. Lineages B.1.617.2 and AY.1-124 are classified as the variants of concern, and are referred to as Delta in accordance with the WHO classification.

According to a number of reports, mutation N501Y in the S gene, found in the genomes of tree discussed above VOC out of four, affects the S protein-ACE2 binding affinity, which may result in increased transmissibility [14–16]. The presence of mutation E484K can allow the virus escape from monoclonal antibodies produced in response to the previous infection with SARS-CoV-2 containing no such mutation, or in response to vaccination [17–19].

It was shown for individuals infected with the Alpha lineage, that the mean duration of the acute infection was longer compared to the values obtained for the strains of other lineages (13.3 days vs. 8.2 days, respectively) [20]. To date, there is no evidence that the SARS-CoV-2 Alpha genetic variant has the potential to escape from monoclonal antibodies, produced in response to vaccination [8, 21–23]. However, according to literature, the South African and Brazilian variants could reduce the efficacy of some vaccines against COVID-19 [7, 24]. Structural analysis of the RBD mutation L452R and the P681R mutation, located in the furin cleavage site, has shown that these mutations may help to increase the efficacy of SARS-CoV-2 binding to the ACE2 receptor and the S1-S2 cleavage rate, thus contributing to the increased transmissibility. The presence of this mutation in RBD indicates the reduced binding with the selected monoclonal antibodies, and may affect their neutralization potency [13, 25].

Based on the foregoing, detection of various coronavirus variants and tracing the spread of the SARS-CoV-2 isolates, being the variants of utmost epidemiological significance, is an important challenge.

Currently, viral genome sequencing is the most common technique, allowing one to discover the rapidly emerging SARS-CoV-2 variants. However, studies involving the use of genome sequencing are expensive and time-consuming. Moreover, the method is also limited by the more stringent requirements to viral load in the sample compared to PCR. Reverse transcription polymerase chain reaction (RT-PCR) is a method, which enables screening the samples containing the SARS-CoV-2 RNA for the already known functionally significant mutations in the S gene of coronavirus, identification of which allows one to define and differentiate the viral VOC of utmost epidemiological significance. The use of oligonucleotide probes with modified LNA nucleotides enables the successful differentiation of even single nucleotide substitutions in the assayed gene sequence [26]. RT-PCR is recommended to increase the output and performance of screening for the top-priority genetic variants by the guidelines on the detection and identification of SARS-CoV-2 variants, issued by the experts of the WHO and ECDC [27]. There are reagent kits of foreign production for detection of important mutations in the gene encoding the S protein by RT-PCR. For example, reagent kit, manufactured by the Korean company PowerChek, identifies mutations, typical for the lineages Delta, Alpha, Beta, and Gamma, and the kit, produced by the French company ID, detects mutations L452R, E484K, and E484Q. In Russia only the reagent kits for detection of mutations typical for the Alpha lineage have been registered (manufactured by the Central Research Institute of Epidemiology and DNA-Technology LLC), the other lineages are not covered.

Thus, the development of the domestic reagent kit for the timely monitoring of the top-priority SARS-CoV-2 variants spread is extremely relevant. For this purpose we have developed the RT-PCR-based method and the reagent kit for detection of the range of mutations in the S gene of coronavirus, typical for four genetic variants of utmost epidemiological significance, as well as for detection of RNA of the SARS-CoV-2 variants Delta, Alpha, and Beta/Gamma (without differentiation between the latter) by identification of the relevant combination of key mutations. The study was aimed to develop the reagent kit allowing one to identify mutations N501Y, P681H and the deletion 69-70del, the combination of which is typical for the SARS-CoV-2 Alpha genetic lineage, mutation E484K, which in combination with the mutation N501Y is typical for the Beta and Gamma lineages, mutations P681R and L452R, the

**Table.** The RNA targets with corresponding fluorescent probes to be detected for each RT-PCR mixture of the kit

Dye of fluorescent probe	RNA target to be detected		
	FAM	JOE	ROX
2-A	mutation N501Y in the S gene of SARS-CoV-2	fragment of codon 501 of the S gene of SARS-CoV-2 with no mutation N501Y	deletion 69-70del in the S gene of SARS-CoV-2
2-B	fragment of codon 484 of the S gene of SARS-CoV-2 with no mutation E484K	mutation E484K in the S gene of SARS-CoV-2	mutation E484Q in the S gene of SARS-CoV-2
2-C	mutations P681H/R in the S gene of SARS-CoV-2	mutation L452R in the S gene of SARS-CoV-2	mutation P681R in the S gene of SARS-CoV-2

combination of which is typical for the Delta lineage, as well as to detect the presence of each of the listed above mutations. The kit may be used as a monitoring tool enabling the prompt acquisition of data on the dynamic changes in the spread of the top-priority SARS-CoV-2 strains (or genetic variants), carrying the functionally significant mutations.

## METHODS

### Design of oligonucleotides for amplification of the S gene target regions

The whole genome alignment of the SARS-CoV-2 sequences, deposited in the GISAID database (accessed 05.10.2021), was used to select the diagnostic primers and probes for detection of mutations, typical for the most epidemiologically significant SARS-CoV-2 genetic lineages. Based on the alignment data, regions carrying the relevant mutations were defined in the viral S gene, primers and LNA-modified probes were selected to detect the following mutations in the genome of the SARS-CoV-2 coronavirus: N501Y, deletions 69-70del, P681H and P681R, L452R, E484K, E484Q. Oligonucleotides were selected in accordance with the standard requirements for selection of primers and TaqMan probes [28, 29] with the use of the Oligo Calc [30] and OligoAnalyzer Tool [31] online resources. Thermodynamic properties and secondary structure of fluorescent probes were assessed using the Mfold Web Server [32]. Oligonucleotides were synthesized by the Genterra (Russia).

### Clinical samples

Nasopharyngeal and oropharyngeal swabs (n = 10,297) were obtained at the Head Center for Hygiene and Epidemiology of the Federal Medical Biological Agency from patients with symptoms of ARVI from December 2020 to September 2021. Swab samples were collected in the test tubes, containing transport medium (Central Research Institute of Epidemiology; Russia), and stored at a temperature of -70 °C before testing.

### RNA extraction

RNA was extracted using the AmpliTest RIBO-prep kit (FSBI SPC FMBA; Russia) in accordance with the manufacturer's guidelines. RNA was eluted in 50 µl of RNA elution buffer. RNA samples were stored at a temperature not exceeding -70 °C before using.

### RT-PCR

Multiplex RT-PCR was carried out in accordance with the instructions. In addition to RNA samples, each test included the positive control samples (PCS) A, B, C, PCS W for the wild type, and the negative control.

Each sample was tested for the presence of mutations, typical for the VOC Alpha, Delta, Beta/Gamma, in three test tubes, containing mixtures 2-A, 2-B, and 2-C. Fluorescence signal accumulation curves for three different channels were analyzed (Table). The volume of RNA solution to be tested was 10 µl.

Amplification program included the following thermal cycling steps: 50 °C — 15 min, 95 °C — 15 min. The following steps were repeated in 45 cycles: 95 °C — 15 s, 60 °C — 30 s, 72 °C — 15 s.

The AmpliTest SARS-CoV-2 VOC v.2 reagent kit was optimized using the Rotor-Gene Q (Qiagen; Germany), CFX96 (Bio-Rad Laboratories; USA), and DTprime (DNA-Technology; Russia) systems.


### Positive control samples (PCS)

The mixture of recombinant plasmids containing the amplified target S-gene fragments, being the targets of the selected primers and probes, was used as a PCS.

The PCR products were purified using the MiniElute Gel Extraction Kit (Qiagen; Germany), ligated into the pGEM-T plasmid vector (Promega; USA), and transformed in *Escherichia coli*. Recombinant plasmids of individual clones were purified using the Plasmid Miniprep Kit (Axygen; USA). The PCS nucleotide sequences were confirmed by Sanger sequencing (cycle sequencing) using the ABI PRISM Big Dye v.3.1 kit (Thermo Fisher Scientific; USA) in accordance with the manufacturer's guidelines, and the Applied Biosystems Sanger Sequencing 3500 Series Genetic Analyzer (Thermo Fisher Scientific; USA). The forward and reverse PCR primers, flanking the fragment to be amplified, were used for sequencing.

The sequencing reaction was carried out in 5 µL of reaction mixture in the 0.2 mL thin wall microtubes. The reaction mixture contained 0.8 µL of specific primer (forward or reverse) with a concentration of 1 µM and 1 µL of the Ready Reaction BigDye Terminator v3.1 mixture (Thermo Fisher Scientific; USA). The components were mixed carefully, and thermal cycling was performed in the SimpliAmp VeriFlex 96 thermal cycler (Thermo Fisher Scientific; USA) in accordance with the following program: initial denaturation 2 min at 95 °C, then 45 cycles (denaturation 20 s at 95 °C, primer annealing 30 s at 60 °C, elongation 1 min at 68 °C). Capillary electrophoresis was performed in the Applied Biosystems Sanger Sequencing 3500 Genetic Analyzer (Thermo Fisher Scientific; USA) after the sample purification of excess dideoxynucleotides and salts.

Four PCSs were developed for the reagent kit: PCS A, containing the fragment of the SARS-CoV-2 S-gene codon 501 carrying the mutation N501Y, and the S-gene fragment carrying the deletion 69-70del; PCS B, containing the mixture of recombinant plasmids with the fragments of the SARS-CoV-2 S-gene codon 484 carrying mutations E484K and E484Q; PCS C, containing plasmids carrying mutations P681H, L452R

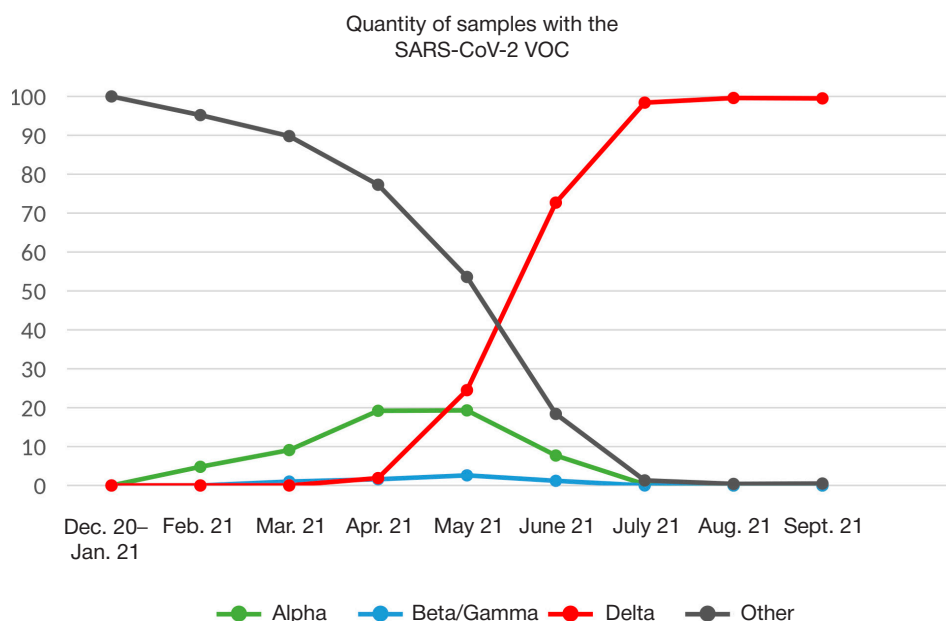

**Fig.** The results for the SARS-CoV-2 VOC spread monitoring performed from December 2020 to September 2021. The samples (total  $n = 10,297$ ) were collected in Moscow and Moscow Region

and P681R in the SARS-CoV-2 S gene; PCS W, containing plasmids with the fragments of the SARS-CoV-2 S-gene codons 484 and 501 with no mutations N501Y and E484K.

The concentration of each PCS was measured by Droplet Digital PCR (ddPCR) with the use of the QX200 Droplet Digital PCR System (Bio-Rad Laboratories; USA). PCSs were introduced at the stage of RT-PCR as the separate samples.

#### Analytical sensitivity and specificity

To assess the analytical sensitivity, we used the SARS-CoV-2 RNA samples, obtained by extraction from biomaterial (smears from the mucous membrane of nasopharynx and oropharynx), which were studied using the diagnostic system developed in order to identify the minimum dilution, at which the samples were detected as positive. The concentration of each sample was previously measured by ddPCR using the QX200 Droplet Digital PCR System (Bio-Rad Laboratories; USA). The sensitivity threshold was defined based on the minimum dilution detected in three iterations.

To assess the analytical specificity of the reagent kit, we used RNA of the SARS-CoV-2 strain No. GK2020/1 from the collection of the N. F. Gamaleya National Research Center, as well as RNA of the strains of coronaviruses HCoV 229E, HCoV OC43, HCoV NI63, SARS-CoV HKU39849, MERS-CoV (European Virus Archive Global 011N-03868 — Coronavirus RNA specificity panel), strains of influenza A virus (H1N1) (ATCC VR-1469), influenza A virus (H3N2) (ATCC VR-776) and influenza B virus (Victoria Lineage) (ATCC VR-1930) from the American Type Culture Collection (ATCC; USA), strains of *Streptococcus pneumoniae* (№ 131116), *Streptococcus pyogenes* (№ 130001), *Haemophilus influenzae* (№ 151221), *Staphylococcus aureus* (№ 201108), *Klebsiella pneumoniae* from the collection of pathogenic microorganisms of the Scientific Centre for Expert Evaluation of Medicinal Products of the Ministry of Health of the Russian Federation, strains of human parainfluenza viruses types 1, 2, 3, human rhinoviruses types 13, 17, 26 from the State Collection of the Smorodintsev Research Institute of Influenza in a concentration not less than  $1 \times 10^6$  GE/mL.

#### Whole genome sequencing (WGS)

Whole genome sequencing (WGS) was used as a reference method for identification of mutations in the RNA samples.

Samples with the cycle threshold (Ct) in RT-PCR not exceeding 25 were used for sequencing. Reverse transcription of the RNA samples was performed using the Reverta-L reagent kit (Centre for Strategic Planning of FMBA; Russia). The resulting cDNA was amplified using the AmpliSeq for Illumina SARS-CoV-2 Research Panel (Illumina; USA), containing 247 amplicons in two pools, which covered the entire SARS-CoV-2 genome. Libraries for the high-throughput sequencing were prepared using the AmpliSeq Library PLUS kit (Illumina; USA). The quality of libraries was evaluated by capillary electrophoresis with the Agilent Bioanalyzer 2100 (Agilent; USA). The library concentrations were measured with the Qubit 4 Fluorometer (Thermo Fisher Scientific; USA) using the HS Qubit dsDNA reagent kit (Thermo Fisher Scientific; USA). Sequencing was carried out with the use of the Illumina NextSeq 550 system (Illumina; USA) and the NextSeq 500/550 v2.5 kit (300 cycles) (Illumina; USA). All the procedures were carried out in accordance with the guidelines, issued by the manufacturers.

#### Statistical analysis

The diagnostic sensitivity (percentage of positive test results in the group of positive samples) was calculated as  $A/(A+B) \times 100\%$ , and the diagnostic specificity (percentage of negative test results in the group of negative samples) was calculated as  $A/(A+B) \times 100\%$ , where  $B$  was the number of assayed samples, which tested positive using the kit, of the total number of samples, and  $A$  was the number of assayed samples, which tested negative using the kit, of the total number of samples with the true positive results. The 95% confidence intervals for the diagnostic characteristics were calculated by the Clopper-Pearson method.

#### RESULTS

When developing the reagent kit, we have selected the SARS-CoV-2 S-gene regions, flanking mutations, typical for the



SARS-CoV-2 genetic lineages Alpha, Beta/Gamma, and Delta, to be the targets. We have selected primers and probes for identification of mutations N501Y, P681H and the deletion 69–70del, the combination of which is typical for the SARS-CoV-2 Alpha genetic lineage (B.1.1.7), mutations P681R and L452R, the combination of which is typical for the Delta genetic lineage (B.1.617.2 and all the AY variants), mutation E484K, which in combination with N501Y is typical for the Beta and Gamma genetic lineages (B.1.351/P.1, these variants have not been differentiated).

The reagent kit has been developed in the format of using three RT-PCR mixtures for analysis. Each mixture enables the detection of three markers, involving registration in three different fluorescence channels. In addition to the key mutations in the S gene, the range of genetic markers, identified using the kit, includes two codons, E484 and N501, matching with the reference SARS-CoV-2 RNA sequence. These serve as endogenous controls for the exclusion of false positive results.

The samples of all the listed above viruses and bacteria were tested in order to define the analytical specificity, no false positive results were obtained. When assaying RNA of the SARS-CoV-2 strain No. GK2020/1 from the collection of the N. F. Gamaleya National Research Center using the kit, no mutations were found in the codons 501, 484, 681, and 69–70 of the SARS-CoV-2 S gene, the results corresponded to the whole genome nucleotide sequence of this isolate, deposited in the GISAID database (EPI\_ISL\_421275). The following result was obtained for all samples containing no SARS-CoV-2 RNA: the samples have no SARS-CoV-2 RNA for analysis.

To assess the analytical sensitivity, we used the biomaterial samples with the known concentration of the SARS-CoV-2 viral RNA, defined by ddPCR. The genetic variants were reproducibly detected up to the SARS-CoV-2 RNA concentration of  $1 \times 10^3$  copies/mL, which defined the analytical sensitivity of the reagent kit.

The diagnostic characteristics of the reagent kit developed were assessed in the clinical trial, involving studying the SARS-CoV-2 RNA samples with genotypes defined by WGS or Sanger sequencing of the S gene. To assess the diagnostic sensitivity, a total of 192 nasopharyngeal and oropharyngeal swab samples were analyzed, including 32 samples containing the SARS-CoV-2 RNA with a combination of mutations in the S gene, characteristic of the Alpha genetic lineage (B.1.1.7), 28 samples containing the SARS-CoV-2 RNA with a combination of mutations, characteristic of the Beta/Gamma genetic lineages (B.1.351/P1), 50 samples containing the SARS-CoV-2 RNA of the Delta genetic lineage (B.1.617.2 and all the AY variants). To assess the diagnostic specificity, a total of 82 RNA samples obtained from nasopharyngeal and oropharyngeal swabs were analyzed, which contained the SARS-CoV-2 RNA with no mutations under analysis. All the target mutations and their combinations, corresponding to the Delta, Alpha or Beta/Gamma genetic lineages, were identified in all samples. No discordant (nonspecific) results were registered. Thus, the diagnostic sensitivity (DS) of 100% was defined (DS with the 95% level of confidence was 94.2–100% for the Delta lineage, 91.1–100% for the Alpha lineage, 89.9–100% for the Beta/Gamma lineages). The diagnostic specificity was 100% (with 95% level of confidence it was 96.4–100%).

The further practical use of the reagent kit showed full compliance of VOC and their key mutations' identification results with the results of WGS performed in 1500 samples of viral RNA of the genetic lineages Alpha, Beta, Gamma, and others (collected from March to May 2021), as well as Delta (collected from June to August 2021). All the samples carried

the relevant combinations of mutations in the S gene, based on which the following SARS-CoV-2 genetic lineages were defined: Delta ( $n = 750$ ), Alpha ( $n = 302$ ), Beta/Gamma ( $n = 32$ ), and others ( $n = 416$ ). The RT-PCR results were completely matched with the data of WGS, so high diagnostic performance of the reagent kit was shown.

From December 1, 2020 to September 30, 2021 more than 10,000 clinical samples of the SARS-CoV-2 RNA were tested, collected in Moscow and Moscow Region within the framework of monitoring the spread of the SARS-CoV-2 VOC. The findings were split per month and presented as a percentage of the variants identified (Fig.). During the studied period (December 2020 to September 2021) we observed the emergence of VOC or single functionally significant mutations and their subsequent replacement by the novel SARS-CoV-2 virus variants. Thus, in the samples, obtained from December 2020 to January 2021, no isolates were found carrying mutations, typical for VOC. During the period the gradual increase in the proportion of distinct lineages was observed. In five months the proportion of the SARS-CoV-2 Beta/Gamma lineages samples changed from 0 to 1.6%, and the proportion of samples belonging to the Alpha lineage increased from 0 to 19.2%. The proportion of the Delta variant rose rapidly to 25% by the end of May, reached 87% by mid-June 2021, exceeded 95% in the first week of July, and was above 99% during the next months until the end of the studied period.

## DISCUSSION

The emergence of the novel SARS-CoV-2 variants remains a grave worldwide concern, since the new variants can have the potential of higher transmissibility, affect the disease duration and severity, reduce the vaccine efficacy, and increase the mortality rate [6, 7, 14–20, 24].

Currently, WGS is the most widely used versatile method for identification of the new SARS-CoV-2 variants. Unfortunately, the method is very time-consuming and expensive. RT-PCR is an efficient alternative method, suitable for identification of the previously defined mutations, the markers of the key lineages, the detection of which is important for monitoring the spread of the genetic variants, and understanding the epidemiological situation.

The outcome of the study is the real time RT-PCR diagnostic test system for identification of the functionally significant mutations allowing one to define the SARS-CoV-2 VOC. The analysis includes the following steps: RNA extraction, RNA reverse transcription, cDNA amplification with the real time fluorescence hybridization detection, and data interpretation.

The kit provided allows one to distinguish RNA of the SARS-CoV-2 Alpha, Beta/Gamma and Delta variants from the other genetic lineages based on the combinations of typical mutations. According to a number of studies [15–23], the presence of distinct mutations (E484K, N501Y, P681H, P681R), associated with different genetic variants of coronavirus, results in the functionally significant alterations in the S-protein structure, including those contributing to the increased transmissibility, reduced neutralization efficiency of antibodies, produced in response to infection with the earlier circulating variants of the virus or in response to vaccination. That is why the detection of individual mutations, without taking into account the combinations of mutations, is also of some epidemiological significance.

The kit shows high analytical sensitivity ( $1 \times 10^3$  copies/mL of the SARS-CoV-2 RNA) for each of the genetic variants to be detected, and 100% analytical specificity in the tested panel of

microorganisms. The diagnostic specificity with the 95% level of confidence is 94.2–100% for the Delta lineage, 91.1–100% for the Alpha lineage, 89.9–100% for the Beta/Gamma lineages, and the diagnostic specificity is 96.4–100% with the 95% level of confidence.

The use of the reagent kit developed within the framework of epidemiological monitoring made it possible to detect the emergence of the SARS-CoV-2 Delta variant in Moscow in April 2021 in a timely manner, and to report promptly the dramatic increase in the proportion of Delta strains among the SARS-CoV-2 variants, detected in the surveyed patients in Moscow Region. It is important to note that the emergence and further spread of the Delta variant in Russia have resulted in the almost complete substitution of the other coronavirus lineages.

## References

1. Bogoch II, Watts A, Thomas-Bachli A, Huber C, Kraemer MUG, Khan K, et al. Pneumonia of unknown aetiology in Wuhan, China: potential for international spread via commercial air travel. *J Travel Med.* 2020; 27 (2): 1–3.
2. Hui DS, I Azhar E, Madani TA, Ntoumi F, Kock R, Dar O, et al. The continuing 2019-nCoV epidemic threat of novel coronaviruses to global health — The latest 2019 novel coronavirus outbreak in Wuhan, China. *Int J Infect Dis.* 2020; 91: 264–6.
3. Rothan HA, Byrareddy SN. The epidemiology and pathogenesis of coronavirus disease (COVID-19) outbreak. *J Autoimmun.* 2020 May; 109:102433.
4. Sánchez CM, Gebauer F, Suñé C, Mendez A, Dopazo J, Enjuanes L. Genetic evolution and tropism of transmissible gastroenteritis coronaviruses. *Virology.* 1992 Sep; 190 (1): 92–105.
5. GISAID EpiFlu™ Database. Available from: <http://www.GISAID.org>.
6. World Health Organization. Available from: <https://www.who.int>.
7. Luring AS, Hodcroft EB. Genetic variants of SARS-CoV-2 — What do they mean? *JAMA.* 2021; 325: 529–31.
8. CDC 2021. Emerging SARS-CoV-2 variants. Available from: <https://www.cdc.gov/coronavirus/2019-ncov/more/science-and-research/scientific-brief-emerging-variants.html>.
9. Grubaugh ND, Hodcroft EB, Fauver JR, Phelan AL, Cevik M. Public health actions to control new SARS-CoV-2 variants. *Cell.* 2021;184 (5): 1127–32.
10. Tegally H, Wilkinson E, Giovanetti M, Iranzadeh A, Fonseca V, Giandhari J, et al. Emergence and rapid spread of a new severe acute respiratory syndrome-related coronavirus 2 (SARS-CoV-2) lineage with multiple spike mutations in South Africa. *MedRxiv.* 2020; 12.21.20248640.
11. Nelson G, Buzko O, Spilman PR, Niazi K, Rabizadeh S, Soon-Shiong PR. Molecular dynamic simulation reveals E484K mutation enhances spike RBD-ACE2 affinity and the combination of E484K, K417N and N501Y mutations (501Y.V2 variant) induces conformational change greater than N501Y mutant alone, potentially resulting in an escape mutant. *bioRxiv.* 2021.01.13.426558.
12. National Institute of Infectious Diseases of JAPAN Brief report 12.01.2021. Available from: <https://www.niid.go.jp/niid/images/epi/corona/covid19-33-en-210112.pdf>.
13. Cherian S, Potdar V, Jadhav S, Yadav P, Gupta N, Das M, et al. Convergent evolution of SARS-CoV-2 spike mutations, L452R, E484Q and P681R, in the second wave of COVID-19 in Maharashtra, India. *Microorganisms.* 2021; 9 (7): 1542.
14. Lacobucci G. Covid-19: New UK variant may be linked to increased death rate, early data indicate. *BMJ.* 2021; 372: n230.
15. Santos JC, Passos GA. The high infectivity of SARS-CoV-2 B.1.1.7 is associated with increased interaction force between Spike-ACE2 caused by the viral N501Y mutation. *bioRxiv.* 2020.12.29.424708.
16. Davies NG, Abbott S, Barnard RC, Jarvis CI, Kucharski AJ, Munday JD, et al. Estimated transmissibility and impact of SARS-CoV-2 lineage B.1.1.7 in England. *Science.* 2021; 372 (6538): eabg3055.
17. Cele S, Gazy I, Jackson L, Hwa SH, Tegally H, Lustig G, et al. Escape of SARS-CoV-2 501Y.V2 from neutralization by convalescent plasma. *Nature.* 2021; 593 (7857): 142–6.
18. Wibmer CK, Ayres F, Hermanus T, Madzivhandila M, Kgagudi P, Oosthuysen B, et al. SARS-CoV-2 501Y.V2 escapes neutralization by South African COVID-19 donor plasma. *Nat Med.* 2021; 27: 622–5.
19. Garcia-Beltran WF, Lam EC, St Denis K, Nitido AD, Garcia ZH, Hauser BM, et al. Multiple SARS-CoV-2 variants escape neutralization by vaccine-induced humoral immunity. *Cell.* 2021; 184 (9): 2372–2383.e9.
20. Kissler S, Fauver JR, Mack C, Tai C, Breban MI, et al. Densely sampled viral trajectories suggest longer duration of acute infection with B.1.1.7 variant relative to non-B.1.1.7 SARS-CoV-2. Harvard University's DASH repository 2021. *medRxiv.* 2021.
21. Callaway E, Mallapaty S. Novavax offers first evidence that COVID vaccines protect people against variants. *Nature.* 2021; 590: 17.
22. Greaney AJ, Loes AN, Crawford KHD, Starr TN, Malone KD, Chu HY, et al. Comprehensive mapping of mutations to the SARS-CoV-2 receptor-binding domain that affect recognition by polyclonal human serum antibodies. *Cell Host Microbe.* 2021; 29 (3): 463–476.e6.
23. Wu K, Werner AP, Moliva JL, Koch M, Choi A, et al. mRNA-1273 vaccine induces neutralizing antibodies against spike mutants from global SARS-CoV-2 variants. *BioRxiv.* 2021. 01.25.427948.
24. Elbe S, Buckland-Merrett G. Data, disease and diplomacy: GISAID's innovative contribution to global health. *Global challenges (Hoboken, NJ).* 2017; 1: 33–46.
25. Starr TN, Greaney AJ, Hilton SK, Ellis D, Crawford KHD, Dingens AS, et al. Deep mutational scanning of SARS-CoV-2 receptor binding domain reveals constraints on folding and ACE2 binding. *Cell.* 2020; 182 (5): 1295–1310.e20.
26. You Y, Moreira BG, Behlke MA, Owczarzy R. Design of LNA probes that improve mismatch discrimination. *Nucleic Acids Res.* 2006; 34 (8): e60.
27. European Centre for Disease Prevention and Control, WHO Regional Office for Europe. Methods for the detection and identification of SARS-CoV-2 variants. 3 March 2021. ECDC and WHO Regional Office for Europe: Stockholm and Copenhagen. 2021.
28. Van Pelt-Verkuil E, van Belkum A, Hays JP. Principles and Technical Aspects of PCR Amplification. Springer Science & Business Media. 2008; 330 p.
29. Yuryev A. Methods in Molecular Biology: PCR Primer Design. Humana Press. 2007; 432 p.
30. Kibbe WA. OligoCalc: an online oligonucleotide properties calculator. *Nucleic Acids Res.* 2007. Available from: <http://biotools.nubic.northwestern.edu/OligoCalc.html>.
31. Integrated DNA Technologies. OligoAnalyzer Tool. Available from:

## CONCLUSIONS

The highly sensitive, specific, and easy to use AmpliTest SARS-CoV-2 VOC v.2 reagent kit for identification of mutations in the S gene of coronavirus, typical for the Alpha, Beta/Gamma and Delta genetic lineages, by RT-PCR has been developed in the Centre for Strategic Planning of the Federal Medical Biological Agency. The kit was validated using the samples containing the coronavirus RNA with genotypes, determined by the SARS-CoV-2 genome sequencing. The PCR results showed a perfect match with the sequencing data. The use of the PCR reagent kit enables fast and efficient assessment of the spread of the coronavirus genetic variants, and immediately taking the appropriate anti-epidemic measures based on the data obtained.

<https://www.idtdna.com/pages/tools/oligoanalyzer>.

32. The mfold Web Server (Hosted by The RNA Institute, College

of Arts and Sciences). Available from: <http://unafold.rna.albany.edu/?q=mfold/DNA-Folding-Form>.

## Литература

1. Bogoch II, Watts A, Thomas-Bachli A, Huber C, Kraemer MUG, Khan K, et al. Pneumonia of unknown aetiology in Wuhan, China: potential for international spread via commercial air travel. *J Travel Med.* 2020; 27 (2): 1–3.
2. Hui DS, I Azhar E, Madani TA, Ntoumi F, Kock R, Dar O, et al. The continuing 2019-nCoV epidemic threat of novel coronaviruses to global health — The latest 2019 novel coronavirus outbreak in Wuhan, China. *Int J Infect Dis.* 2020; 91: 264–6.
3. Rothan HA, Byrareddy SN. The epidemiology and pathogenesis of coronavirus disease (COVID-19) outbreak. *J Autoimmun.* 2020 May; 109:102433.
4. Sánchez CM, Gebauer F, Suñé C, Mendez A, Dopazo J, Enjuanes L. Genetic evolution and tropism of transmissible gastroenteritis coronaviruses. *Virology.* 1992 Sep; 190 (1): 92–105.
5. GISAID EpiFlu™ Database. Available from: <http://www.GISAID.org>.
6. World Health Organization. Available from: <https://www.who.int>.
7. Luring AS, Hodcroft EB. Genetic variants of SARS-CoV-2 — What do they mean? *JAMA.* 2021; 325: 529–31.
8. CDC 2021. Emerging SARS-CoV-2 variants. Available from: <https://www.cdc.gov/coronavirus/2019-ncov/more/science-and-research/scientific-brief-emerging-variants.html>.
9. Grubaugh ND, Hodcroft EB, Fauver JR, Phelan AL, Cevik M. Public health actions to control new SARS-CoV-2 variants. *Cell.* 2021;184 (5): 1127–32.
10. Tegally H, Wilkinson E, Giovanetti M, Iranzadeh A, Fonseca V, Giandhari J, et al. Emergence and rapid spread of a new severe acute respiratory syndrome-related coronavirus 2 (SARS-CoV-2) lineage with multiple spike mutations in South Africa. *MedRxiv.* 2020; 12.21.20248640.
11. Nelson G, Buzko O, Spilman PR, Niazi K, Rabizadeh S, Soon-Shiong PR. Molecular dynamic simulation reveals E484K mutation enhances spike RBD-ACE2 affinity and the combination of E484K, K417N and N501Y mutations (501Y.V2 variant) induces conformational change greater than N501Y mutant alone, potentially resulting in an escape mutant. *bioRxiv.* 2021.01.13.426558.
12. National Institute of Infectious Diseases of JAPAN Brief report 12.01.2021. Available from: <https://www.niid.go.jp/niid/images/epi/corona/covid19-33-en-210112.pdf>.
13. Cherian S, Potdar V, Jadhav S, Yadav P, Gupta N, Das M, et al. Convergent evolution of SARS-CoV-2 spike mutations, L452R, E484Q and P681R, in the second wave of COVID-19 in Maharashtra, India. *Microorganisms.* 2021; 9 (7): 1542.
14. Lacobucci G. Covid-19: New UK variant may be linked to increased death rate, early data indicate. *BMJ.* 2021; 372: n230.
15. Santos JC, Passos GA. The high infectivity of SARS-CoV-2 B.1.1.7 is associated with increased interaction force between Spike-ACE2 caused by the viral N501Y mutation. *bioRxiv.* 2020.12.29.424708.
16. Davies NG, Abbott S, Barnard RC, Jarvis CI, Kucharski AJ, Munday JD, et al. Estimated transmissibility and impact of SARS-CoV-2 lineage B.1.1.7 in England. *Science.* 2021; 372 (6538): eabg3055.
17. Cele S, Gazy I, Jackson L, Hwa SH, Tegally H, Lustig G, et al. Escape of SARS-CoV-2 501Y.V2 from neutralization by convalescent plasma. *Nature.* 2021; 593 (7857): 142–6.
18. Wibmer CK, Ayres F, Hermanus T, Madzivhandila M, Kgagudi P, Oosthuysen B, et al. SARS-CoV-2 501Y.V2 escapes neutralization by South African COVID-19 donor plasma. *Nat Med.* 2021; 27: 622–5.
19. Garcia-Beltran WF, Lam EC, St Denis K, Nitido AD, Garcia ZH, Hauser BM, et al. Multiple SARS-CoV-2 variants escape neutralization by vaccine-induced humoral immunity. *Cell.* 2021; 184 (9): 2372–2383.e9.
20. Kissler S, Fauver JR, Mack C, Tai C, Breban MI, et al. Densely sampled viral trajectories suggest longer duration of acute infection with B.1.1.7 variant relative to non-B.1.1.7 SARS-CoV-2. *Harvard University's DASH repository.* 2021.
21. Callaway E, Mallapaty S. Novavax offers first evidence that COVID vaccines protect people against variants. *Nature.* 2021; 590: 17.
22. Greaney AJ, Loes AN, Crawford KHD, Starr TN, Malone KD, Chu HY, et al. Comprehensive mapping of mutations to the SARS-CoV-2 receptor-binding domain that affect recognition by polyclonal human serum antibodies. *Cell Host Microbe.* 2021; 29 (3): 463–476.e6.
23. Wu K, Werner AP, Moliva JI, Koch M, Choi A, et al. mRNA-1273 vaccine induces neutralizing antibodies against spike mutants from global SARS-CoV-2 variants. *BioRxiv.* 2021. 01.25.427948.
24. Elbe S, Buckland-Merrett G. Data, disease and diplomacy: GISAID's innovative contribution to global health. *Global challenges (Hoboken, NJ).* 2017; 1: 33–46.
25. Starr TN, Greaney AJ, Hilton SK, Ellis D, Crawford KHD, Dingens AS, et al. Deep mutational scanning of SARS-CoV-2 receptor binding domain reveals constraints on folding and ACE2 binding. *Cell.* 2020; 182 (5): 1295–1310.e20.
26. You Y, Moreira BG, Behlke MA, Owczarzy R. Design of LNA probes that improve mismatch discrimination. *Nucleic Acids Res.* 2006; 34 (8): e60.
27. European Centre for Disease Prevention and Control, WHO Regional Office for Europe. Methods for the detection and identification of SARS-CoV-2 variants. 3 March 2021. ECDC and WHO Regional Office for Europe: Stockholm and Copenhagen. 2021.
28. Van Pelt-Verkuil E, van Belkum A, Hays JP. Principles and Technical Aspects of PCR Amplification. Springer Science & Business Media. 2008; 330 p.
29. Yuryev A. Methods in Molecular Biology: PCR Primer Design. Humana Press. 2007; 432 p.
30. Kibbe WA. OligoCalc: an online oligonucleotide properties calculator. *Nucleic Acids Res.* 2007. Available from: <http://biotools.nubic.northwestern.edu/OligoCalc.html>.
31. Integrated DNA Technologies. OligoAnalyzer Tool. Available from: <https://www.idtdna.com/pages/tools/oligoanalyzer>.
32. The mfold Web Server (Hosted by The RNA Institute, College of Arts and Sciences). Available from: <http://unafold.rna.albany.edu/?q=mfold/DNA-Folding-Form>.